



## Elephant-II FANS Flexible, Suction Ureteral Access Sheath


The Elephant-II Flexible and Navigable Suction (FANS) ureteral access sheath has a flexible 10cm distal portion, which is directed into the desired calyx by the flexible ureteroscope. Stone dust and fragments are evacuated via suction connected to the side arm, controlled by a thumb switch.



**Atraumatic:** The distal tip is designed with a 3mm section that consists solely of soft PTFE and Pebax®, excluding any metal coil, to protect tissue from damage when advancing and navigating the sheath.



**Integrated Self-Seal:** The innovative, conical seal guides the scope into the sheath without damaging the deflecting section and can be released quickly should a large basketed stone need to be removed.

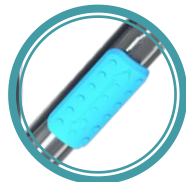

**Slider Suction:** Easily control the suction level with a convenient thumb slider, leaving your hands free to control the ureteroscope and sheath.



**Scope Savings:** Purchase the Elephant II sheath at a discounted rate together with our range of HugeMed single use flexible ureteroscopes.

Size / Length	45cm ♀	55cm ♂
9/11 Fr	YGE-AGS09045	YGE-AGS09055
10/12 Fr	YGE-AGS10045	YGE-AGS10055
11/13 Fr	YGE-AGS11045	YGE-AGS11055
12/14 Fr	YGE-AGS12035	YGE-AGS12055

