



**JENA SURGICAL**


**novus med**

SURGICAL INNOVATION, KNOWLEDGE & RELIABILITY

# JENA Surgical 150w MultiPulse HoPLUS Ho:YAG Laser



**Optimised for HoLEP:** The high power 150w Ho:YAG laser was the first on the market to reach a 100Hz repetition rate with 0.2 - 6 J pulse energy. The EPS (Enhanced Power Setting) and long pulse duration of 1700  $\mu$ s ensures excellent haemostasis during HoLEP.

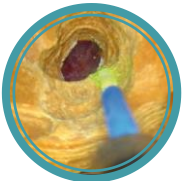

**Efficient Stone Dusting:** The high repetition rate of 100 Hz and long pulse duration of up to 1700  $\mu$ s not only facilitates efficient stone dusting but also decreases the retropulsion effect, resulting in quicker procedures.



**32-Amp Power Supply:** Unlike other 100+ watt lasers that require a 64-amp power supply, the MultiPulse HoPLUS generates 150 watts of power from a single phase 32-amp power supply, which is available in most UK theatres.




**Made in Germany:** Based in the Jena Optical Valley in Germany, JenaSurgical, boasts over 40 years of expertise in the development and distribution of laser systems for multi-disciplinary surgical applications.



**Available via NHS Supply Chain:** Available on the Operating Theatres Equipment and Related Accessories and Services (Reference: Project\_972), Lot 4 – Medical Lasers and Related Accessories and Services.

Product Code	Description
J19-1150-000 819-1154-100	MultiPulse HoPLUS 150w / 100Hz
J19-1150-000	MultiPulse HoPLUS 100w / 100Hz
J19-1181-000 819-1154-075	MultiPulse HoPLUS 100w / 75Hz
J19-1181-000	MultiPulse HoPLUS 80w / 75Hz

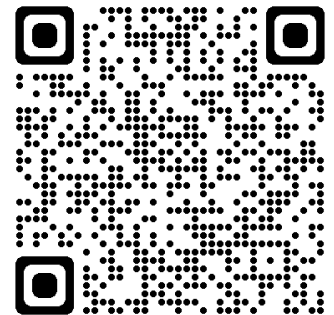

**novus med**  
SURGICAL INNOVATION, KNOWLEDGE & RELIABILITY

+44 (0) 1926 895089  
customerservice@novus-med.com

