

HD

HugeMed

novus med

SURGICAL INNOVATION, KNOWLEDGE & RELIABILITY



Single-use Range of Flexible Cystoscopes



Portable: Ideal for difficult catheterisations, the 8” HD display monitor can be easily moved between departments and is mains powered or has a 4hr battery time. It features 32GB memory and USB connectivity for external storage.



Versatile: The 8” monitor can be connected to an external monitor via HDMI cable or for a more permanent solution, our HUV-01 control unit can be integrated into an existing camera system in theatres or outpatients.



High Definition: The latest generation of CMOS sensor offers HD visualisation. Whilst the bullet-tip design leads to a smoother insertion and increased patient comfort.




Suction as Standard: The suction function helps clear the visual field by removing fluid and debris providing better visualisation of the area being examined during cystoscopy.



Choice of Sizes: Available in the standard 15Fr version with a 6.6Fr channel or a slimmer 12Fr scope for tighter anatomy / PCNL, with a 3.6Fr channel. A Tuohy Borst adaptor is included on the instrument port as standard.

Product Code	Description
CY-M52H	10x 15Fr HD flexible cystoscope, 6.6Fr working channel
CY-M40H	10x 12Fr HD flexible cystoscope, 3.6Fr working channel
VLM-03	HD 8” Portable 1024x768 pixels display monitor
HUV-01	HD Control unit, 1024 x 768 pixels resolution

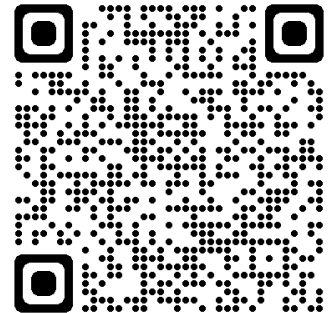

novus med
SURGICAL INNOVATION, KNOWLEDGE & RELIABILITY

LG1.SUFC.V1

+44 (0) 1926 895089
customerservice@novus-med.com

