



ELEPHANT-II Flexible, Suction Ureteral Access Sheath


The ELEPHANT-II Flexible and Navigable Suction Ureteral Access Sheath has a flexible 10cm distal portion, which is directed into the desired calyx by the flexible ureteroscope. Stone dust and fragments are evacuated via suction connected to the side arm, controlled by a thumb switch.



Atraumatic: The distal tip is designed with a 3mm section that consists solely of soft PTFE and Pebax®, excluding any metal coil, to protect tissue from damage when advancing and navigating the sheath.



Integrated Self-Seal: The innovative, conical seal guides the scope into the sheath without damaging the deflecting section and can be released quickly should a large basketed stone need to be removed.

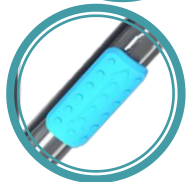

Slider Suction: Easily control the suction level with a convenient thumb slider, leaving your hands free to control the ureteroscope and sheath.



Stone Trap Included: A new, smaller stone collection pot reduces plastic waste and is more ergonomic when connected to the sheath.

Available on



Supply Chain

UoM	Size / Length	40cm	50cm
Box of 5	9/11 Fr	YGE-AGS09045 NPC: FUX85237	YGE-AGS09055 NPC: FUX85244
	10/12 Fr	YGE-AGS10045 NPC: FUX85243	YGE-AGS10055 NPC: FUX85242
	11/13 Fr	YGE-AGS11045 NPC: FUX85241	YGE-AGS11055 NPC: FUX85240
	12/14 Fr	YGE-AGS12045 NPC: FUX85238	YGE-AGS12055 NPC: FUX85239